Journal of Food and Drug Analysis, Vol. 16, No. 3, 2008, Pages 39-45

# Determination of Aldicarb, Propoxur, Carbofuran, Carbaryl and Methiocarb Residues in Honey by HPLC with Post-column Derivatization and Fluorescence Detection after Elution from a Florisil Column

FERIDE KOC<sup>1\*</sup>, YUSUF YIGIT<sup>2</sup>, YAVUZ KURSAD DAS<sup>3</sup>, YASEMIN GUREL<sup>2</sup> AND CEVDET YARALI<sup>4</sup>

<sup>1.</sup> Department of Pharmacology and Toxicology, Faculty of Veterinary Medicine, Atatürk University, Erzurum 25400, Turkey

<sup>2.</sup> Central Veterinary Control and Research Institute, Toxicology Laboratory, Ankara, Turkey

<sup>3</sup> Department of Pharmacology and Toxicology, Faculty of Veterinary Medicine, 19 Mayis University, Samsun, Turkey

<sup>4.</sup> Central Veterinary Control and Research Institute Genetic Laboratory, Ankara, Turkey

(Received: march 20, 2007; Accepted: September 16, 2007)

### ABSTRACT

The aim of this study was to develop an easy analytical method for determining aldicarb, propoxur, carbofuran, carbaryl and methiocarb residues in honey by high performance liquid chromatography with post-column derivatization and fluorescence detection (HPLC/PCD-FL) for routine analysis. A florisil packed column was used in sample preparation step and showed good performance. The recoveries ranged between 72.02 and 92.02% for various spiking levels (50, 100, 200 ng/g). The limits of detection were 4–5 ng/g, depending on different *N*-methylcarbamate pesticides. The extraction method used in the present study was easy, inexpensive, and rapid. Thus, it could be used efficiently in the routine monitoring studies.

Key words: honey, liquid chromatography, post-column derivatization, N-methylcarbamate pesticide, residue, florisil

## **INTRODUCTION**

Carbamate pesticides, a class of highly effective commercial pesticides, have been used worldwide since the beginning of 1960s due to their high insecticide and nematocide effects. Their structures are N-substituted carbamic acid esters<sup>(1)</sup>. Carbamate pesticides began to replace organochlorine such as DDT and organophosphorous pesticides, due to their low environmental persistence and low toxic effects on mammalians $^{(2,3)}$ . However, carbamate pesticides are suspected carcinogens and mutagens<sup>(4)</sup> because of they are also acetylcholinesterase inhibitors. The increasing use of carbamate pesticides poses a risk to apiculture and human environment. Thus, it is necessary to quantify their residue amounts in foods to prevent harmful effects on animals, human and the environment<sup>(4-</sup> <sup>9)</sup>. Honey, a natural product of bees, is composed primarily of sugars and water. Quality of honey is directly related to its floral origin and the region of production. Its quality control is performed to monitor residues<sup>(10)</sup>. According to the Directive 96/23/EC, it is necessary to research N-meth-

\* Author for correspondence. Tel: +90-442-631-4193 ext. 1007; Fax: +90-442-631-4188; E-mail: feridekoc@yahoo.com ylcarbamate pesticide residues in honey<sup>(11)</sup>. Many methods are available for determination of carbamate pesticide residues, such as ELISA, cholinesterase-based biosensors and chromatography<sup>(12,13)</sup>. For chromatography, different sample preparation methods have been applied to the determination of carbamate pesticides, including liquidliquid extraction (LLE), solid-phase extraction (SPE), and supercritical fluid extraction (SFE). SPE requires much lower volume in organic solvent usage<sup>(14-20)</sup>. On the other hand, common techniques for determining carbamate pesticides are based on chromatographic techniques, such as gas chromatography (GC), high-performance liquid chromatography (HPLC), liquid chromatography mass spectrometer (LC-MS), GC-MS, LC-MS-MS and GC-MS-MS<sup>(21-25)</sup>. For polar, less volatile and thermally labile pesticides, such as carbamate pesticides, HPLC is obviously the prefered approach because of its high sensitivity and selectivity<sup>(25-31)</sup>. Derivatization reaction is employed mainly to increase sensitivity or selectivity for the detector, particularly fluorescence detector. Among the precolumn and post-column derivatization methods reviewed, post-column method is the most commonly used<sup>(32-39)</sup>.

The carbamate hydrolysis and derivatization reac-

### tion is shown in Figure 1.

In this study, we preferred to use HPLC post-column derivatization and fluorescence detection and aimed to develop a new, simple multiresidue method for determining several *N*-methylcarbamate pesticide residues in honey by using a florisil packed-column. The florisil was used in previous studies<sup>(40-44)</sup>. The method described in this study was rapid and suitable for screening pesticide residues in monitoring programs. Clean-up was performed using activated florisil packed-column as SPE technique. No methods have so far been documented for *N*-methylcarbamate pesticides in honey by HPLC with post-column derivatization and fluorescence detection after elution from activated florisil.


# MATERIALS AND METHODS


#### I. Biolocigal Materials

Honey samples were collected from local markets in Turkey and used to develop method for *N*-methylcarbamate pesticide residue analyses. Samples were stored during work at -20°C.

#### II. Chemicals, Reagents and Standards

In this study, all solvents (acetonitrile, methanol, n-hexane) were of analytical grade and obtained from Merck (Darmstadt, Germany). Hydrolysis reagent, thiofluor, [o-phthalaldehyde (OPA)] was supplied from Pickering Laboratories (California, USA). The standards of *N*-methylcarbamate pesticide (aldicarb 98.5%, propoxur 99.0% carbofuran 98.0%, carbaryl 99.5% and methiocarb 98.5% purity) were purchased from Dr. Ehrenstorfer GmbH (Augsburg, Germany). Florisil (60-100 mesh) and dichloromethane were obtained from J.T. Baker (Deventer, Netherlands). Anhydrous sodium sulfate was from Merck (Darmstadt, Germany). Ultra pure water was supplied by Elga (London, UK).

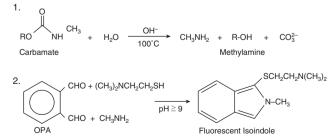

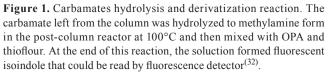

Journal of Food and Drug Analysis, Vol. 16, No. 3, 2008

# III. Instruments

HPLC/PCD-FL was used for analyses. The HPLC system consisted of an Agilent 1100 Series, including gradient pump, fluorescence detection system, autosampler, vacuum membrane degasser, and software Chemstation Rev.b01.03 (2001-2005) (Avondale, PA, USA). The post column system was supplied by Prometheus (Athena, Greece).

Chromatographic separation was achieved on a reversed-phase C18 column, 5  $\mu$ m, 150  $\times$  4.6 mm I.D obtained from Pickering Laboratories (California, USA). The guard column XDB-C<sub>8</sub> (4.6  $\times$  12.5 mm I.D.) was obtained from Agilent (Avondale, PA, USA).

#### IV. Chromatographic Conditions

The mobile phase consisted of a mixture of water and acetonitrile at a ratio of 90:10 (v/v). The time program of mobile phase system during analysis period is shown in Table 1.

The injection volume was 20  $\mu$ L. Monitoring wavelength was set at exication 330 nm, and emission 465 nm. Flow-rate of mobile phase was set at 1 mL/min, and the temperature of the column was 42°C. Hydrolysis reaction temperature was set at 100°C and the flow-rate of hydrolysis and OPA reagents were both 0.3 mL/min<sup>(45)</sup>.

# (I) Preparation of the Reagents

The hydrolysis reagent (950 mL) was prepared in No. 1 bottle.

The OPA diluent (950 mL) was prepared in No. 2 bottle.

## (II) Preparation of the OPA Reagent

Five milliliter of the OPA diluent was taken from No. 2 bottle and 2 g of thiofluor was dissolved in this OPA diluent. One hundred milligram of OPA was dissolved in 10 mL of HPLC-grade methanol. Both of solutions were added to No. 2 bottle.

#### V. Standards and Working Solutions

Concentration of stock standard solution was 10 ng/

Table 1. Time program of mobile phase system

| Step | Time (min) | Interval | Water, % | ACN, %* |
|------|------------|----------|----------|---------|
| 0    | 0          | 0        | 90       | 10      |
| 1    | 0-35       | 35       | 90       | 10      |
| 2    | 35-40      | 5        | 10       | 90      |
| 3    | 40-60      | 20       | 90       | 10      |

\*ACN: Acetonitrile.

Journal of Food and Drug Analysis, Vol. 16, No. 3, 2008

 $\mu$ L. Working standard solutions of all carbamate pesticides were prepared at six different concentrations (0.01, 0.05, 0.1, 0.25, 0.5, 1 ng/ $\mu$ L) in methanol. Calibration curve was prepared in the range of 0.01–1 ng/ $\mu$ L for all carbamate pesticides (aldicarb, propoxur, carbofuran, carbaryl, methiocarb) (n = 6).

# VI. Sample Preparation

Frozen honey samples were thawed at room temperature (25°C) and honey (1 g) sample was weighed in a 25 mL glass tube and added with 2 mL of methanol. The sample was mixed by handling. Florisil was activated at 130°C for at least 48 hr. A little cotton wool was put into a glass column (200 × 20 mm I.D.) Then, 12 g prepared florisil was added with anhydrous sodium sulfate (10 g). After mixing with 20 mL of *n*-hexane/dichloromethane (1/1), the mixture was transferred to the glass column. The column was conditioned. The sample was applied to the glass column and eluted with 15 mL of n-hexane/ dichloromethane mixture (1/1, v/v). The eluant was collected in a glass tube and then evaporated (at 40°C, reduced pressure). Residue was dissolved with 1 mL of acetonitrile.

### VII. Validation

The following parameters were determined for validation of developed method: linearity, precision, accuracy, limit of detection (LOD), limit of quantitation (LOQ), and recovery.

For linearity, six levels  $(0.01-1 \text{ ng/}\mu\text{L})$  of calibration series with six injections at each concentration level were determined. Calibration curve was calculated automatically using software.


Recoveries were calculated as the percentage using extraction process (with activated florisil) after spiking from 50 to 200 ng/g with three different levels of carbamate pesticides (n = 6).

Precision and accuracy were defined based on spiked pesticide-free honey samples at three concentrations by calibration curve prepared daily (n = 6). Relative standard deviation (RSD) and relative standard error (accuracy) were determined. All samples were freshly prepared from reference standard solutions (10 ng/µL). Accuracy was calculated according to this formula: (Concentration found – concentration known) × 100 / Concentration known.

Considering a signal-to-noise ratio of 3 and 10 for LOD and LOQ, respectively.

## **RESULTS AND DISCUSSION**

Chromatogram of carbamate pesticide mix standard (1 ng/ $\mu$ L) is shown in Figure 2. Retention times were 16.528, 18.942, 19.758, 20.437 and 23.972 min for aldi-

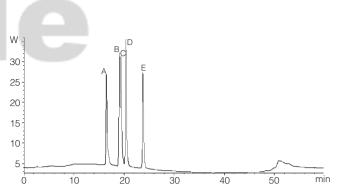

**Figure 2.** HPLC/PCD-FL chromatogram on C18 column by injecting 20  $\mu$ L of standard solution (1 ng/ $\mu$ L): A: aldicarb peak, B: Propoxur peak, C: Carbofuran peak, D: Carbaryl peak, E: Methiocarb peak.

**Table 2.** Limit of detections (LOD, ng/g), limit of quantitations (LOQ, ng/g) and correlation cofficients ( $R^2$ ) in 0.01-1 ng/ $\mu$ L standard solution range with six levels

| Pesticide  | LOD (ng/g) | LOQ (ng/g) | Correlation (R <sup>2</sup> ) |
|------------|------------|------------|-------------------------------|
| Aldicarb   | 5          | 15         | 0.9998                        |
| Propoxur   | 4          | 16         | 0.9999                        |
| Carbofuran | 4          | 16         | 0.9999                        |
| Carbaryl   | 4          | 16         | 0.9999                        |
| Methiocarb | 4          | 16         | 0.9999                        |

carb, propoxur, carbofuran, carbaryl and methiocarb, respectively (Figure 2).

Calibration graphs were established within the range of 0.01–1 ng/ $\mu$ L with correlation of coefficients from 0.9998 to 0.9999 for all analytes. The correlation of coefficients (R<sup>2</sup>), LODs and LOQs are shown in Table 2. LOD of aldicarb was 5 ng/g whereas LODs of others were 4 ng/g. LOQs were found nearly 3.72 times of LODs.

Blasco *et al.*<sup>(46)</sup> reported that Porteguese honey was more contaminated than Spanish honey. However, little pesticide residues were found in Porteguese honey.

Khan *et al.*<sup>(47)</sup> reported that among 36 honey samples collected from different bee colonies, carbaryl was detected in one honey sample each of Apis mellifera (901.8 ppb) and Apis florae (800.0 ppb) in India.

Miyata *et al.*<sup>(48)</sup> reported that the extract was purified on a silica gel column and GC was used for determination of carbamate pesticides. No pesticide residue was found in the 56 samples of polished rice from Thailand. The detection limit of the proposed method was 0.01 ppm for carbamate pesticides. In the present study, LODs were 4-5 ng/g. We used HPLC/PCD/FL because it was more sensitive and selective than GC for carbamate pesticides. Due to their physical chemical properties such as thermal instability and polarity, carbamate pesticides are difficult or even impossible to be analysed using GC techniques without the time-consuming derivatization process. Due to lower sensitivity of LC/UV, HPLC with fluorescence detection by post-column derivatization is 42

the most widely used method for the analysis of carbamate pesticides in  $foods^{(49,50)}$ . LODs obtained from the present study were consistent with other papers<sup>(51,52)</sup>.

For recovery studies, honey samples were spiked by adding 50, 100, and 200 ng/g standards of *N*-methylcarbamate pesticide. Chromatogram of blank sample and spiked samples (200 ng/g) on C18 column are shown in Figures 3 and 4, respectively.

Precision (RSD) and accuracy ranged 1.77-9.23% and -7.98-27.98%, respectively (n = 6).

Recoveries of the spiked samples, RSD and accuracy are shown in Table 3.

Average recoveries of the spiked samples were 80.09, 77.34, 85.11, 78.02, 79.78% for aldicarb, propoxur, carbofuran, carbaryl and methiocarb, respectively. The recoveries were between 72.02 and 92.02% for varione spiking levels (50, 100, 200 ng/g). Bernal *et al.*<sup>(53)</sup> and Kawasaki *et al.*<sup>(54)</sup> reported that recovery of SPE florisil packed-column was higher than liquid-liquid extraction method. Recoveries in the present study were consistent with previous studies<sup>(55-57)</sup>.

Previous studies reported that the extraction of pesticide residues from honey (SPE with florisil packed column) was easier and cheaper than the commonly used procedures<sup>(42,44,57)</sup>. Jimenez *et al.*<sup>(42,58)</sup> reported that florisil cartridges contributes to reducing the matrix effects. In this study, matrix effect was not observed.


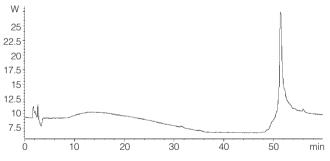


Figure 3. Chromatogram of blank honey sample on C18 column.

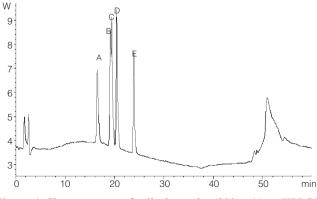

**Figure 4.** Chromatogram of spiked samples (200 ng/g) on HPLC/ PCD–FL. A: aldicarb peak, B: Propoxur peak, C: Carbofuran peak, D: Carbaryl peak, E: Methiocarb peak.

| Pesticide  | Spiked level<br>(ng/g)<br>( n = 6) | Recovery<br>(%) | RSD<br>(precision, %) | Relative error<br>(accuracy, % ) |
|------------|------------------------------------|-----------------|-----------------------|----------------------------------|
| Aldicarb   | 50                                 | 82.11           | 5.20                  | -17.89                           |
|            | 100                                | 83.28           | 7.36                  | -16.72                           |
|            | 200                                | 74.89           | 4.22                  | -25.11                           |
|            | average                            | 80.09           | 8.90                  | -19.91                           |
| Propoxur   | 50                                 | 74.94           | 1.92                  | -25.06                           |
|            | 100                                | 80.08           | 9.23                  | -19.92                           |
|            | 200                                | 77.01           | 4.39                  | -22.99                           |
|            | average                            | 77.34           | 6.42                  | -22.66                           |
| Carbofuran | 50                                 | 78.22           | 2.76                  | -21.78                           |
|            | 100                                | 92.02           | 9.11                  | -07.98                           |
|            | 200                                | 85.10           | 4.95                  | -14.90                           |
|            | average                            | 85.11           | 7.28                  | -14.89                           |
| Carbaryl   | 50                                 | 77.09           | 1.77                  | -22.91                           |
|            | 100                                | 76.52           | 8.23                  | -23.48                           |
|            | 200                                | 80.47           | 5.19                  | -19.53                           |
|            | average                            | 78.02           | 7.94                  | -21.98                           |
| Methiocarb | 50                                 | 72.02           | 3.51                  | -27.98                           |
|            | 100                                | 86.66           | 8.40                  | -13.34                           |
|            | 200                                | 80.68           | 6.68                  | -19.32                           |
|            | Average                            | 79.78           | 7.55                  | -20.22                           |

Journal of Food and Drug Analysis, Vol. 16, No. 3, 2008

**Table 3.** Recoveries of spiked samples (n = 6), relative standard deviations (RSD, %) and accuracy in 3 levels

This situation resulted from activated florisil and anhydrous sodium sulfate.

# CONCLUSIONS

In the present study, a new method was developed by HPLC with post-column derivatization and fluorescence detection after elution from a florisil column for several carbamate pesticides. A simultaneous extraction/clean-up method was proposed for determination of aldicarb, propoxur, carbofuran, carbaryl and methiocarb residues in honey using a mixture of florisil and anhydrous sodium sulfate in small glass columns. The florisil column is capable of affecting isolation and clean up of these pesticide residues from honey samples. The main advantage of the developed method is that extraction and cleanup are performed in one step with small volume of organic solvent. Additionally, the validation parameters show that the developed method is a sensitive, reproducible and reliable alternative to the normally used methods. Moreover this inexpensive, easy and rapid method

Journal of Food and Drug Analysis, Vol. 16, No. 3, 2008

can be used in the routine and monitoring studies. The analysis can be confirmed by LC/MS in future studies and the developed method can be used to analyze other carbamate pesticides and their main metabolites.

# REFERENCES

- 1. Soriana, J. M., Jimenez, B., Font, G. and Molto, J. C. 2001. Metabolites in water by solid phase extraction and liquid chromatography: a review. Crit. Rev. Anal. Chem. 31: 19-52.
- Fernandez, J. M., Vazquez, P. P and Vidal, J. L. M. 2000. Analysis of *N*-methylcarbamate pesticides insecticides and some of their main metabolites in urine with liquid chromatography using diode array detection and electrospray mass spectrometry. Anal. Chim. Acta 412: 131-139.
- 3. Petropoulou, S. S. E., Gikas, E., Tsarbopoulos, A. and Siskos, P. A. 2006. Gas chromatographic-tandem mass spectrometric method for the quantitation of carbofuran, carbaryl and their main metabolites in applicators' urine. J. Chromatogr. A 1108: 99-110.
- Zhang, J. and Lee, H. K. 2006. Application of liquidphase microextraction and on-column derivatization combined with gas chromatography-mass spectrometry to the determination of carbamate pesticides. J. Chromatogr. A 1117: 31-37.
- 5. Abad, A., Moreno, M. J., Pelegri, R., Martinez, M. I., Saez, A., Gamon, M. and Montaya, A. 1999. Determination of carbaryl, carbofuran and methiocarb in cucumbers and strawberies by monoclonal enzyme immunoassays and high performance liquid chromatography with fluorescence detection. An analytical comparison. J. Chromatogr. A 833: 3-12.
- Rawn, D. F. K., Roscoe, V., Krakalovich, T. and Hanson, C., 2006. *N*-methyl carbamate concentrations and dietary intake estimates for apple and grape juices available on the retail market in Canada. Food Addit. Contam. 21: 553-563.
- Rawn, D. F. K., Roscoe, V., Trelka, Y., Hanson, C., Krakalovich, T. and Dabeka, R. W. 2006. *N*-methyl carbamate pesticide residues in conventional and organic infant foods available on the Canadian retail market, 2001-03. Food Addit. Contam. 23: 651-659.
- Cabras, P., Tuberoso, C., Melis. M. and Martini, G. 1992. Multiresidue method for pesticide determination in wine by high-performance liquid chromatography. J. Agric. Food Chem. 40: 817-819.
- Inthavong, D. and Bordet, F. 2002. Multiresidue analytical method for carbamate pesticides in meat products analysis by LC/DAD and by LC/ACPI/MS/MS. 4th E. P. R. W. May 28-31, Roma, Italy.
- Luo, X. B., Chen, B. and Yao, S. Z. 2004. High-performance liquid chromatography with electrospray mass spectrometry for rapid and sensitive determination of sanguinarine and chelerythrine in exogenously con-

taminated honey. Chromatographia 60: 347-351.

- ANON. 1996. Council Directive 96/23/EC of 29 April 1996. On measures to monitor certain substances and residues there of in live animals and animal products and repealing Directives 85/358/EEC and 86/469/EEC and Decisions 89/187/EEC and 91/664/EEC. No: L 125/11. 23.05.1996.
- Wang, S., Yu, C. and Wang, J. 2005. Enzyme immunoassay for the determination of carbaryl residues in agricultural products. Food Addit. Contam. 22: 735-742.
- Skladal, P., Nunes, G. S., Yamanaka, H. and Ribeiro, M. L. 1997. Detection of carbamate pesticides in vegetable samples using cholinesterase-based biosensors. Electroanalysis. 9: 1083-1087.
- Sherma, J. and Shafik, T. M. 1975. A multiclass, multiresidue analytical method for determining pesticide residues in air. Arch. Environ. Contam. Toxicol. 3: 55-71.
- 15. Molina, M., Perez-Bendito, D. and Silva, M. 1999. Multi-residue analysis of *N*-methylcarbamate pesticides and their hydrolytic metabolites in environmental waters by use of solid-phase extraction and micellar electrokinetic chromatography. Electrophoresis 20: 3439-3449.
- Santos Delgado, M. J., Rubio, B. S., Fernandez-Tostado, G. T. and Polo-Diez, L. M. 2001. Stability studies of carbamate pesticides and analysis by gas chromatography with flame ionization and nitrogen–phosphorus detection. J. Chromatogr. A 921: 287-296.
- 17. King, J. W. and Zhang, Z. 2002. Derivatization reactions of carbamate pesticides in supercritical carbon dioxide. Anal. Bioanal. Chem. 374: 88-92.
- Carabias-Martinez, R., Garcia-Hermida, U., Rodriguez-Gonzalo, E. and Ruano-Miguel, L. 2005. Behaviour of carbamate pesticides in gas chromatography and their determination with solid-phase extraction and solidphase microextraction as preconcentration steps. J. Sep. Sci. 28: 2130-2138.
- Suma, R., Sarin, R. K., Saiprakash, P. K. and Ramakrishna, S. 2005. Simple liquid chromatographic method for the rapid and simultaneous determination of propoxur and its major metabolite isopropoxy phenol in rat blood and urine using solid-phase extraction. J. Anal. Toxicol. 29: 728-733.
- 20. Pang, G. F., Fan, C. L., Liu, Y. M., Cao, Y. Z., Zhang, J. J., Fu, B. L., Li, X. M., Li, Z.Y. and Wu, Y. P. 2006. Multi-residue method for the determination of 450 pesticide residues in honey, fruit juice and wine by double-cartridge solid-phase extraction/gas chromatog-raphy-mass spectrometry and liquid chromatography-tandem mass spectrometry. Food Addit. Contam. 23: 777-810.
- Pirard, C., Widart, J., Nguyen, B. K., Deleuze, C., Heudt, L. and Haubruge, E. 2007. Development and validation of a multi-residue method for pesticide determination in honey using on-column liquid–liquid extraction and liquid chromatography–tandem mass

spectrometry. J. Chromatogr. A. 1152: 116-123.

- 22. Petropoulou, S. S. E., Tsarbopoulos A. and Siskos, P. A. 2006. determination of carbofuran, carbaryl and their metabolites in plasma samples of agricultural populations using gas chromatography-tandem mass spectrometry. Anal. Biol. Chem. 385: 1444-1456.
- 23. Mayer-Helm, B., Hofbauer, L. and Muller, J. 2006. Development of a multi-residue method for the determination of 18 carbamates in tobacco by high-performance liquid chromatography/positive electrospray ionisation tandem mass spectrometry. Rapid Commun. Mass Spectrom. 20: 529-536.
- Sagratini, G., Manes, J., Giardina, D., Damiani, P. and Yolada, P. 2007. Analysis of carbamate and phenylurea pesticide residues in fruit juices by solid-phase microextraction and liquid chromatography–mass spectrometry. J. Chromatogr. A 1147: 135-143.
- Krause, R. T. 1988. High-performance liquid chromatographic determination of aryl *N*-methylcarbamate residues using post-column hydrolysis electrochemical detection. J. Chromatogr. A 442: 333-343.
- Chiron, S. and Barcelo, D. 1993. Determination of pesticides in drinking water by on-line solid-phase disk extraction followed by various liquid chromatographic systems. J. Chromatogr. 645: 125-134.
- 27. Slobodnik, J., Hoekstra-Oussoren, S. J. F., Jager, M. E., Honing, M., van Baar, B. L. M. and Brinkman, U. A. T. 1996. On-line solid-phase extraction-liquid chromatography-particle beam mass spectrometry and gas chromatography-mass spectrometry of carbamate pesticides. Analyst 121: 1327-1334.
- Möder, M., Popp, P., Eisert, R. and Pawliszyn, J. 1999. Determination of polar pesticides in soil by solid phase microextraction coupled to high-performance liquid chromatography-electrospray/mass spectrometry. Fresenius J. Anal. Chem. 363: 680-685.
- Pico, Y., Font, G., Molto, J. C. and Manes, J. 2000. Pesticide residue determination in fruit and vegetables by liquid chromatography–mass spectrometry. J. Chromatogr. A 882: 153-173.
- 30. Vandecasteele, K., Gaus, I., Debreuck, W. and Walraevens, K. 2000. Identification and quantification of 77 pesticides in groundwater using solid phase coupled to liquid liquid microextraction and reversed-phase liquid chromatography. Anal. Chem. 72: 3093-3101.
- Wu, J., Tragas, C., Lord, H. and Pawliszyn, J. 2002. Analysis of polar pesticides in water and wine samples by automated in-tube solid-phase microextraction coupled with high-performance liquid chromatography-mass spectrometry. J. Chromatogr. A 976: 357-367.
- 32. Application Manual. 1999. Carbamates Version 1.1. Pickering Laboratories, Mountain View. CA. U. S. A.
- Bassett, M. V., Wendelken, S. C., Pepich, B. V. and Munch, D. J. 2003. Improvements to EPA Method 531.1 for the analysis of carbamates that resulted in the development of U.S. EPA Method 531.2. J. Chromatogr. Sci.

Journal of Food and Drug Analysis, Vol. 16, No. 3, 2008

41: 100-106.

- McGarvey, B. D. 1994. Derivatization reactions applicable to pesticide determination by high-performance liquid-chromatography. J. Chromatogr. B-Biomed. Appl. 659: 243-257.
- 35. Tsumura, Y., Ujita, K., Nakamura, Y., Tonogai, Y. and Ito, Y. 1995. Simultaneous determination of aldicarb, ethiofencarb, methiocarb and their oxidized metabolites in grains, fruits and vegetables by high-performance liquid-chromatography. J. Food Protect. 58: 217-222.
- 36. Tsumura, Y., Nakamura, Y., Yoshii, K., Tonogai, Y., Higo, M. and Shibata, T. 1998. Simultaneous determination of aldicarb, ethiofencarb, methiocarb and their oxidized metabolites in grains, fruits and vegetables by high-performance liquid-chromatography. J. Food Hyg. Soc. Jpn. 39: 357-367.
- Hiemstra, M. and De Kok, A. 2007. Comprehensive multi-residue method for the target analysis of pesticides in crops using liquid chromatography-tandem mass spectrometry. J. Chromatogr. A. 1154: 3-25.
- Schenck, F. J., Podhorniak, L. V., Hobbs, J, Casanova, J. and Donoghue, D. 2006. Liquid chromatographic determination of *N*-methyl carbamate pesticide residues at low parts-per-billion levels in eggs. J. AOAC Int. 89: 196-200.
- 39. De Kok, A., Hiemstra, M. and Vreeker, C. P. 1990. Optimization of postcolumn hydrolysis reaction on solid phases for the routine high-performance liquid chromatographic determination of *N*-methylcarbamate pesticides in food products. J. Chromatogr. 507: 459-472.
- 40. Honing, M., Riu, J., Barcelo, D., vanBaar, B. L. M. and Brinkman, U. A. T. 1996. Determination of ten carbamate pesticides in aquatic and sediment samples by liquid chromatography-ionspray and thermospray mass spectrometry. J. Chromatogr. A 733: 283-294.
- 41. Lancas, F. M., Rissato, S. R. And Galhiane, M. S. 1996. Off-line SFE-CZE analysis of carbamate residues in tobacco samples. Chromatographia 42: 323-328.
- 42. Jimenez, J. J., Bernal, J. L., del Nozal., M. J., Toribio, L. and Martin, M. L. 1998. Gas chromatography with electron-capture and nitrogen-phosphorus detection in the analysis of pesticides in honey after elution from a florisil column-influence of the honey matrix on the quantitative results. J. Chromatogr. A 823: 381-387.
- 43. Zhang, W. G., Chu, X. G., Cai, H. X., An, J. and Li, C. J. 2006. Simultaneous determination of 109 pesticides in unpolished rice by a combination of gel permeation chromatography and florisil column purification, and gas chromatography/mass spectrometry. Rapid Commun. Mass Spectrom. 20: 609-617.
- 44. Debertrand, N., Durand, G. and Barcelo, D. 1991. Extraction, cleanup and liquid chromatographic-diode array determination of carbamate pesticides in soil samples. J. Environ. Sci. Heal. A 26: 575-597.
- 45. Pelosi, P., Stefanelli P., Attard Barbini, D., Generali, T., Amendola, G., Girolimetti, S., Vanni, F. and Di

Journal of Food and Drug Analysis, Vol. 16, No. 3, 2008

Muccio, A. 2002. The Italian National Reference Laboratory (Pesticide Residues Section of the ISS-Istituto Superiore di Sanita (National Institute of Health)-Roma, Italy.

- Blasco C., Fernandez, M., Pena, A., Lino, C., Silveira, M. I., Font, G. and Pico, Y. 2003. Assessment of pesticide residues in honey samples from Portugal and Spain. J. Agric. Food Chem. 51: 8132-8138.
- 47. Khan, M. S., Kumari, B., Rohilla, H. F., Kaushik, H. D. and Arora, R. K. 2004. Analysis of insecticide residues in honeys from apiary (Apis mellifera) and wild honey bee (Apis dorsata and apis florae) colonies in India. J. Apicult. Res. 43: 79-82.
- 48. Miyata, M., Kamakura, K., Narita, M., Hirahara, Y., Okamoto, K., Hasegawa, M., Koiguchi, S., Miyoshi, T., Yamana T. and Tonogai, Y. 1994. Studies on systematic determination of organophosphorus, carbamate, organochlorine and pyrethroid pesticides in polished rice by gas-chromatography. J. Food Hyg. Soc. Jpn. 35: 276-286.
- 49. Caballo-Lopez, A. and Luque de Castr, M. D. 2003. Continuous ultrasound-assisted extraction coupled to on line filtration-solid-phase extraction-column liquid chromatography-post column derivatisation-fluorescence detection for the determination of *N*-methylcarbamates in soil and food. J. Chromatogr. A 998: 51-59.
- Liu, M., Hashi, Y., Song, Y. and Lin, J. M. 2005. Simultaneous determination of carbamate and organophosphorus pesticides in fruits and vegetables by liquid chromatography-mass spectrometry. J. Chromatogr. A 1097: 183-187.
- Yu, W., Wang, C. and Chu, X. 1998. Simultaneous determinations of carbamate pesticides in cereal by high performance liquid chromatography with postcolumn fluorescence derivatization. Se. Pu. 16: 430-432.

- 52. Kobayashi, M., Nagayama, T., Takano, I., Tamura, Y., Kimura, N., Kitayama, K. and Saito, K. 2002. Simultaneous determination of *N*-methylcarbamate pesticides and their metabolites in agricultural products. Shokuhin Eiseigaku Zasshi 43: 133-143.
- 53. Bernal, J. L., Jimenez, J. J., del Nozal, M. J., Higes, M. and Llorente, J. 2000. Gas chromatographic determination of acrinathrine and 3-phenoxybenzaldehyde residues in honey. J. Chromatogr. A 882: 239-243.
- Kawasaki, M., Inoue, T., Fukuhara, K. and Uchiyama, S. 1999. Study on GC/MS (SIM) for determination of carbamate and organonitrogen pesticides in foods with simple clean-up by SPE method. J. Food Hyg. Soc. Jpn. 40: 382-390.
- 55. Sanchez-Brunete, C., Albero, B., Miguel, E. and Tadeo, J. L. 2002. Determination of insecticides in honey by matrix solid-phase dispersion and gas chromatography with nitrogen-phosphorus detection and mass spectrometric confirmation. J AOAC Int. 85: 128-133.
- 56. Wang, J., Cheung, W. and Grant, D. 2005. Determination of pesticides in apple-based infant foods using liquid chromatography electrospray ionization tandem mass spectrometry. J. Agric. Food Chem. 53: 528-537.
- Kadenczki, L., Arpad, Z., Gardi, I., Ambrus, A., Gyorfi, L., Reese, G. and Ebing, W. 1992. Column extraction of residues of several pesticides from fruits and vegetables - a simple multiresidue analysis method. J. AOAC Int. 75: 53-61.
- 58. Jimenez, J. J., Bernal, J. L., del Nozal, M. J., Toribio, L. and Arias, E. 2001. Analysis of pesticide residues in wine by solid-phase extraction and gas chromatography with electron capture and nitrogen-phosphorus detection. J. Chromatogr. A 919: 147-156.